

INGLESE - cod.22428/A 05.24



**ORMESA®**  
• MADE IN ITALY SINCE 1980 •

CLIP



# CLIP

## For dynamic families.

CLIP by Ormesa is a **practical and quick seating system, for people with disabilities, children and adults.**

It is **easy to handle, fold and push.**

Its wheels are suitable for all terrain and **folding operations are fast and user friendly.**

Thanks to a **wide range of available components,** it adapts to different needs.

Its **low weight** and **easy folding** make it the ideal solution for those families who are often on the move or who live in cities with required and difficult passages on pavements and crowded areas.



PRODUCT  
MARKED 

CLIP is a **MEDICAL DEVICE** that **COMPLIES WITH THE SAFETY REQUIREMENTS** defined with the **REGULATION (EU) 2017/745**



Clip 1, 2, 3 passed the **CRASH TEST** according to the **ISO 7176/19**

# Lightweight and robust

## Backrest and seat with internal removable aluminum bars



### Double Cross Structure

Located behind the back rest and under the seat in order to guarantee support and comfort.

### Internal bars

**Made in aluminum**, present both in the backrest and seat. They can be shaped to meet specific postural needs.

These bars guarantee the correct stiffness to the back rest and seat, while still maintaining breathability.

**Easily removable for upholstery washing.**



### Back rest

Adjustable in inclination into 3 positions, through zippers. Located behind the backrest, there is a basket for objects.

### Seat

With a slight inclination to increase comfort.



ORMESA



**NEW**

## Adjustable handles

in height. **TOOLS FREE!**



## Removable footrest/s

Extractable, height-adjustable for sizes 1, 2, 3.  
Height-adjustable and foldable for size 4.



**NEW**

## Wheels

puncture-proof, suitable for all types of terrain.

**Clip 1, 2, 3:** front swivel wheels 15 cm;  
20 cm, **with independent brakes.**

**Clip 4:** 17.5 cm front swivel wheels;  
25 cm, **with independent drum brakes.**

**Overcoming step bar.**

# High performance

## Cutting-edge upholstery


CLIP has an **extremely breathable, comfortable, durable, easy-to-clean, and quick-to-dry removable upholstery.**

In the central part of the backrest and seat, the covering is made of **highly innovative and breathable material, both internally and externally.**

The inner fabric is three-dimensional, **lightweight**, and ensures:

- **high air circulation**
- **optimal weight distribution**
- **excellent comfort.**

The outer fabric is also extremely **resistant to abrasion** and is **fire-retardant** according to EN 1021-1:2014 and EN 1021-2:2014 standards.

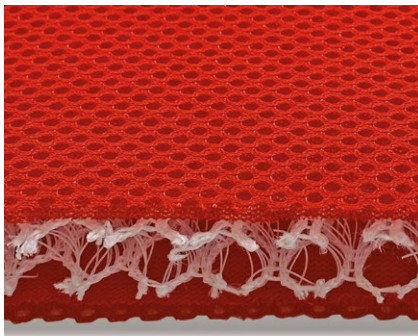

  
**Easily washable  
with fast  
drying times!**



# CLIP 1 • 2 • 3





# CLIP 4



# Components



**835 Headrest with parietal supports.**



**816 Trunk lateral containment pads.**



**853 Vest harness.**



**906 Five point harness.**



**903 Five point vest harness.**



**894 45° pelvic belt.**



**834N Padded and removable abduction block.**

# Components



**858 Shopping basket.**



**827 Foot straps.**



**812 Front wheel directional lock.**



**819 Canopy.**



**825 Rain cover.**  
to fix to the 819 canopy.



**818 Thermal cover.**

### 891 Set of 4 tie down hooks

for the forward-facing transport of the pushchair in a motor vehicle (private cars, buses..). **Only for CLIP 1, 2 and 3.**

See User and Maintenance Manual provided with the aid.






# CLIP 1.2.3



CLIP	A Backrest width	B Overall width	C Overall height	D Backrest inclination	E Seat bottom depth	F Backrest height	G Seat bottom inclination	H ground to seat bottom height	I Seat to footrest height	L Handle to wheel distance	M Overall length	N Folded stroller length	O Folded stroller height	P Folded stroller width	Q User height	Total weight	Max capacity
<b>1</b>	40 cm 15,7 in	51 cm 20 in	101-108 cm 39,7-42,5 in	0° - 12°	25 cm 9,8 in	55 cm 21,6 in	18°	46 cm 18,1 in	12 - 20 cm 4,7 - 7,8 in	78 cm 30,7 in	84 cm 33 in	114 cm 44,8 in	30 cm 11,8 in	34 cm 13,3 in	110 - 125 cm 43,3 - 49,2 in	10,3 kg 22,7 lb	30 kg 66,1 lb
<b>2</b>	40 cm 15,7 in	51 cm 20 in	101-108 cm 39,7-42,5 in	0° - 12°	30 cm 11,8 in	59 cm 23,2 in	18°	48 cm 18,8 in	21 - 31 cm 8,2 - 12,2 in	78 cm 30,7 in	92 cm 36,2 in	114 cm 44,8 in	30 cm 11,8 in	34 cm 13,3 in	120 - 140 cm 47,2 - 55,1 in	10,5 kg 23,1 lb	35 kg 77,1 lb
<b>3</b>	40 cm 15,7 in	51 cm 20 in	106-113 cm 41,7-44,4 in	0° - 12°	35 cm 13,7 in	63 cm 24,8 in	18°	49 cm 19,2 in	25 - 35 cm 9,8 - 13,7 in	84 cm 33 in	98 cm 38,5 in	118 cm 46,4 in	34 cm 13,3 in	36 cm 14 in	135 - 155 cm 53,1 - 61 in	10,7 kg 23,5 lb	40 kg 88,1 lb


# CLIP<sup>4</sup>

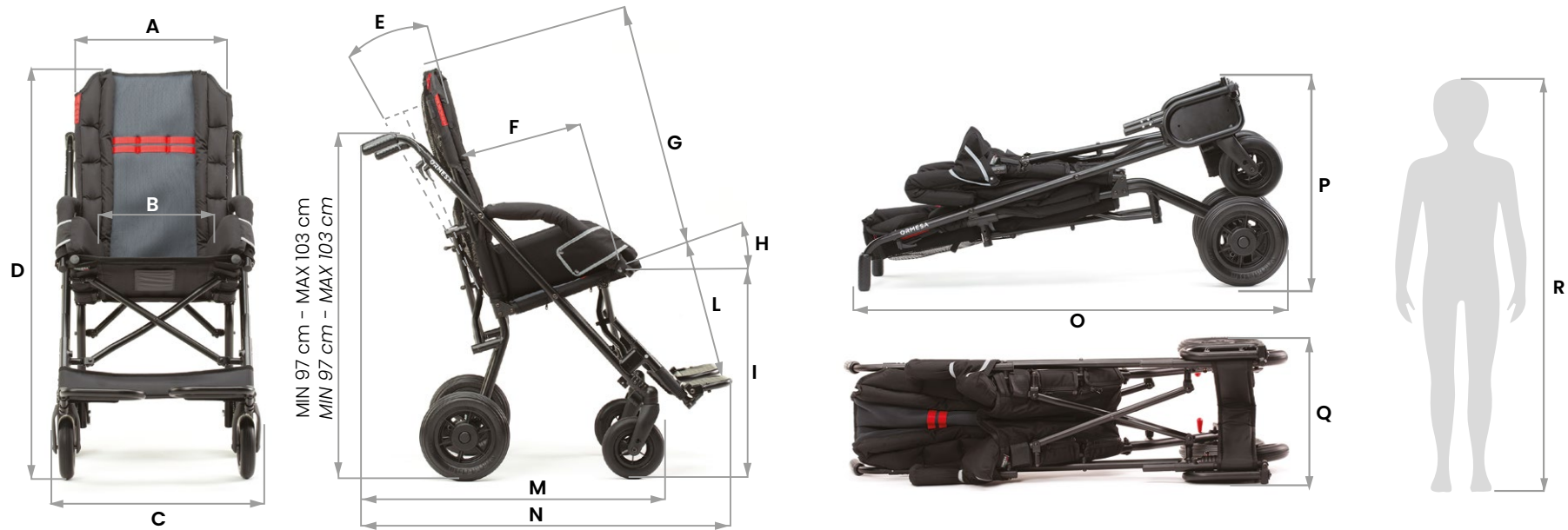

MIN 97 cm - MAX 103 cm  
MIN 97 cm - MAX 103 cm

CLIP	A Backrest width	B Seat bottom width	C Overall width	D Overall height	E Backrest inclination	F Seat bottom depth	G Backrest height	H Seat bottom inclination	I ground to seat bottom height	L Seat to footrest height	M Wheel to handles distance	N Overall length	O Folded stroller length	P Folded stroller height	Q Folded stroller width	R User height	Total weight	Max capacity
4	45 cm 17,7 in	33* cm 12,9* in	57 cm 22,4 in	113 cm 44,4 in	0° - 10°	38 cm 14,9 in	71 cm 27,9 in	20°	59 cm 23,2 in	37 - 47 cm 14,5 - 18,5 in	89 cm 35 in	108 cm 42,5 in	117 cm 46 in	60 cm 23,6 in	37 cm 14,5 in	150 - 170 cm 59 - 66,9 in	13,9 kg 30,6 lb	76 kg 167,5 lb

\*seat width is 45 cm (17,7 in) if you remove the lateral paddings

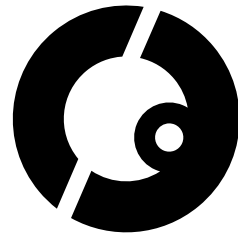

**CE** complies with MDR (EU) 2017/745

COMPANY WITH  
QUALITY SYSTEM  
CERTIFIED BY DNV  
[ISO 13485](#)

---

**Ormesa srl** - Via delle Industrie 6/8 | 06034 Foligno (PG) - ITALY  
Tel. +39 0742 22927 | FAX +39 0742 22637 | [info@ormesa.com](mailto:info@ormesa.com)

**[www.ormesa.com](http://www.ormesa.com)**

