



**ORMESA®**  
• MADE IN ITALY SINCE 1980 •

Duecentotré 5  
**DOCCIA**  
203-5 SHOWER CHAIR



# Duecentotre 5 **DOCCIA**

203-5 SHOWER CHAIR

## Safety in daily hygiene.

**Duecentotre5 Doccia** by Ormesa is a **shower chair for adults**. It is safe and comfortable, designed to let the user be in his semi-lying or sitting position while taking a shower

The **bathroom** is the environment where most **domestic accidents** or **accidental falls** occur: washing a person with motor problems can be uncomfortable and dangerous for both the user and the caregiver. The shower chair **Duecentotre5 Doccia** meets the needs of **safety** and **practicality** of both, ensuring **stability** and **hygiene**.

The **wrapping backrest** has a **"V" shape** for a better containment of the trunk; its structure – together with the **tilted seat** – prevents the user to slide forward and favors the **unloading of the trunk on several points**, excluding the spine from rubbing and compression.

PRODUCT  
MARKED **CE**

203-5 SHOWER CHAIR is a MEDICAL DEVICE COMPLIES WITH THE SAFETY REQUIREMENTS defined with the REGULATION (EU) 2017/745





# Widespread well-being

## Technology for better living

Duecentotre-5 Shower chair has an aluminium frame with a non-toxic epoxy varnish finish.

Its **covering** is made of a special fire-proof, visco-elastic, hygienic, washable and highly wear-resistant **material**. It adapts to the shape of the user's shoulders and glutes and it guarantees a **wide breathability**, as well as an effective water reflux.

The **wrapping backrest** has a **"V" shape** for a better containment of the trunk and it can be fixed in two positions in order to have different depth adjustment of the seat, thanks to the back rack locking.

The **ergonomic seat** has a 15° tilt in order to prevent the User from slipping forward

The **reclining backrest and legrest** allow the caregiver to place the user semi-reclined or seated, keeping the **lateral containment of the trunk-pelvis** constant.



# Technical features and adjustments



## Headrest

with antibacterial coating. It ensures proper head positioning. Adjustable in tilt, height and in forward-backward throw knobs.



## Armrests

extractable and revolving.



## Tip-up footrest

It is adjustable in height and turnable to facilitate the approaching of the care giver and the sitting and standing up of the user.



## Wheels

12 cm diameter. Back wheels with brake.

# Components




**834 Abduction block**  
adjustable in depth.



**894 45° pelvic belt**


# DuecentoTre5 Doccia

## sizes and weight

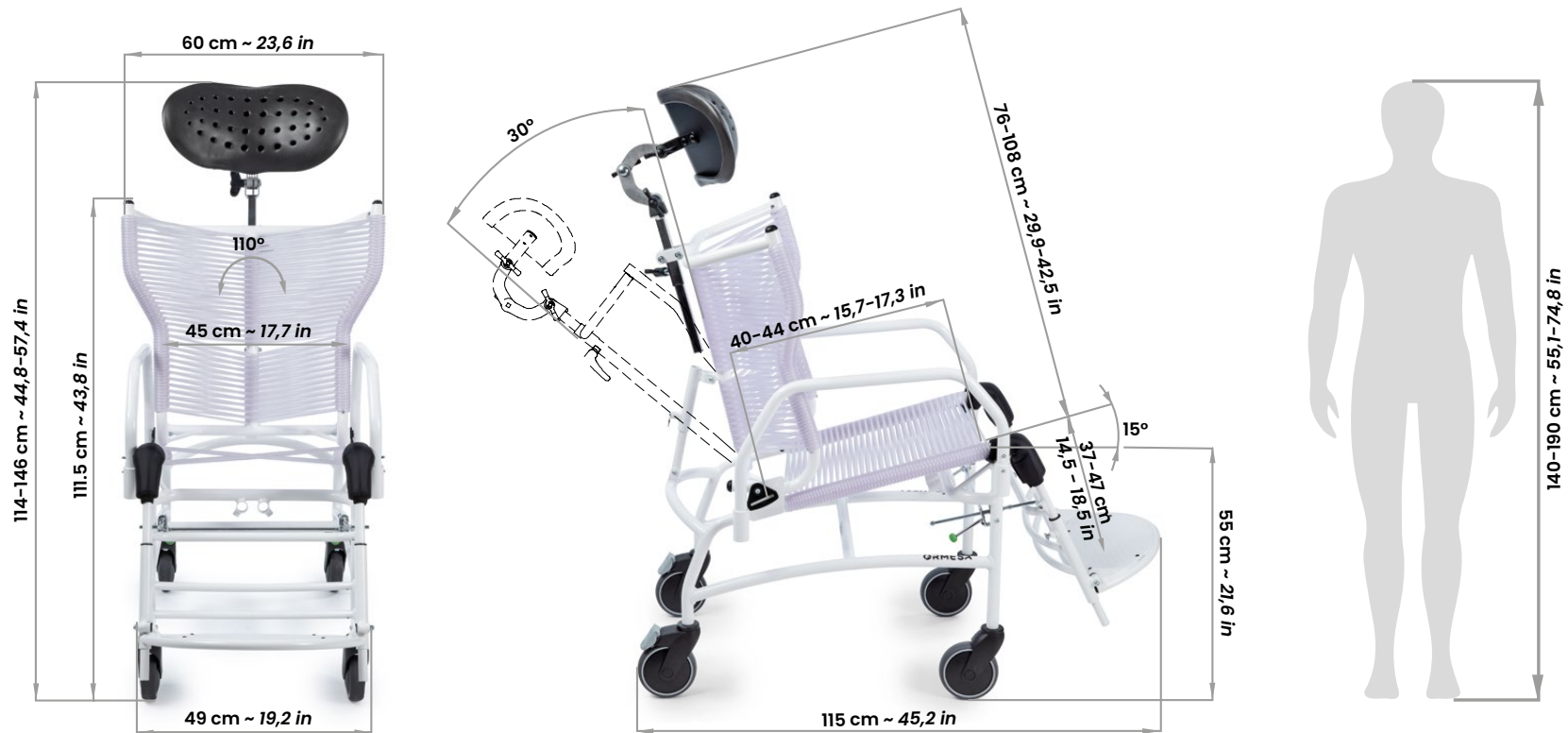

Total Weight 15 Kg ~ 33lb

Max. Load 110 Kg ~ 242,5 lb

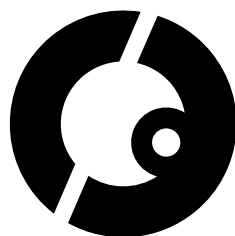

 complies with MDR (EU) 2017/745

COMPANY WITH  
QUALITY SYSTEM  
CERTIFIED BY DNV GL  
= ISO 13485 =

---

**Ormesa srl** - Via delle Industrie 6/8 | 06034 Foligno (PG) - ITALY  
Tel. +39 0742 22927 | FAX +39 0742 22637 | [info@ormesa.com](mailto:info@ormesa.com)

**[www.ormesa.com](http://www.ormesa.com)**