



ORMESA®
• MADE IN ITALY SINCE 1980 •

OBI



OBI

Compact and durable.

OBI by Ormesa is a safe and durable **stroller for children and teens**, equipped with a folding steel frame.

The “Dupont Cordura” nylon seat unit is extremely long-lived and water-washable, **ideal for intensive and continuous use.**

OBI is equipped with a **wide range of components.**



OBI passed **CRASH TEST** according to **ISO 7176/19**

PRODUCT MARKED **CE**

OBI is a **MEDICAL DEVICE** that **COMPLIES WITH THE SAFETY REQUIREMENTS** defined with the **REGULATION (EU) 2017/745**



Technical features and adjustments



Steel folding frame

Painted with non-toxic epoxy varnish.



Removable upholstery

covered with fire retardant fabric, resistant colours and machine washable.



CERTIFICATE
MED annex A.1/3.19-20-21

Nylon seat unit

natural looking "Dupont Cordura", resistant to longtime use and water washable.



Rigid backrest

adjustable in tilt on 4 positions, thanks to a simple quick release mechanism.



Tip-up footrest

with heel pad. It is adjustable in height.



Wheels

Front wheels: 17,5 cm swivel and unpierceable sliding on ball bearings.

Rear wheels: 30 cm unpierceable sliding on ball bearings.



Components



835 Headrest with parietal supports

with two side elements, which can be used one at a time (either the left side or the right side). By turning them upside down, their height position changes.



816 Side supports




834 Abduction block



853 Adjustable harness



894 45° pelvic belt

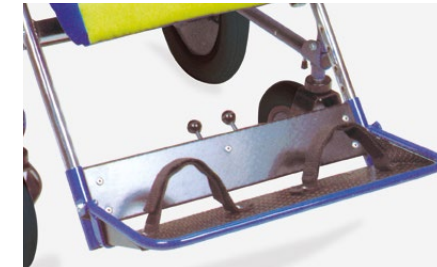

827 Foot straps

Components



839 Front handle



837 Tray without hollow



819 Canopy



825 Rain cover



818 Thermal cover



858 Net basket

Components



812 Directional locks for front wheels




891 Set of 4 tie down hooks for the forward-facing transport of the pushchair in a motor vehicle (private cars, buses..).

See User and Maintenance Manual provided with the aid.



Obi sizes and weight


	obi 2	obi 3	obi 4
A: overall width	60 cm ~ 23,6 in	62,5 cm ~ 24,6 in	66,5 cm ~ 26,1 in
B: seat bottom width	33 cm ~ 12,9 in	37 cm ~ 14,5 in	41 cm ~ 16,1 in
C: overall length	99 cm ~ 38,9 in	109 cm ~ 42,9 in	124 cm ~ 48,8 in
D: overall height	97 cm ~ 38,1 in	103 cm ~ 40,5 in	105 cm ~ 41,3 in
E: seat bottom depth	30 cm ~ 11,8 in	36 cm ~ 14,1 in	40 cm ~ 15,7 in
F: backrest height	62 cm ~ 24,4 in	68 cm ~ 26,7 in	76 cm ~ 29,9 in
G: seat to footrest height	from 23 to 35 cm from 9 to 13,7 in	from 29 to 44 cm from 11,4 to 17,3 in	from 32 to 47 cm from 12,5 to 18,5 in
L: folded stroller length	92 cm ~ 36,2 in	95 cm ~ 37,4 in	103 cm ~ 40,5 in
M: folded stroller height	50 cm ~ 19,6 in	46 cm ~ 18,1 in	51 cm ~ 20 in
H: user height	from 120 to 140 cm from 47,2 to 55,1 in	from 135 to 155 cm from 53,1 to 61 in	from 150 to 175 cm from 59 to 68,8 in
Total weight	17,4 kg ~ 38,3 lb	18,5 kg ~ 40,7 lb	20,1 kg ~ 44,3 lb
Max capacity	50 kg ~ 110,2 lb	60 kg ~ 132,2 lb	70 kg ~ 154,3 lb

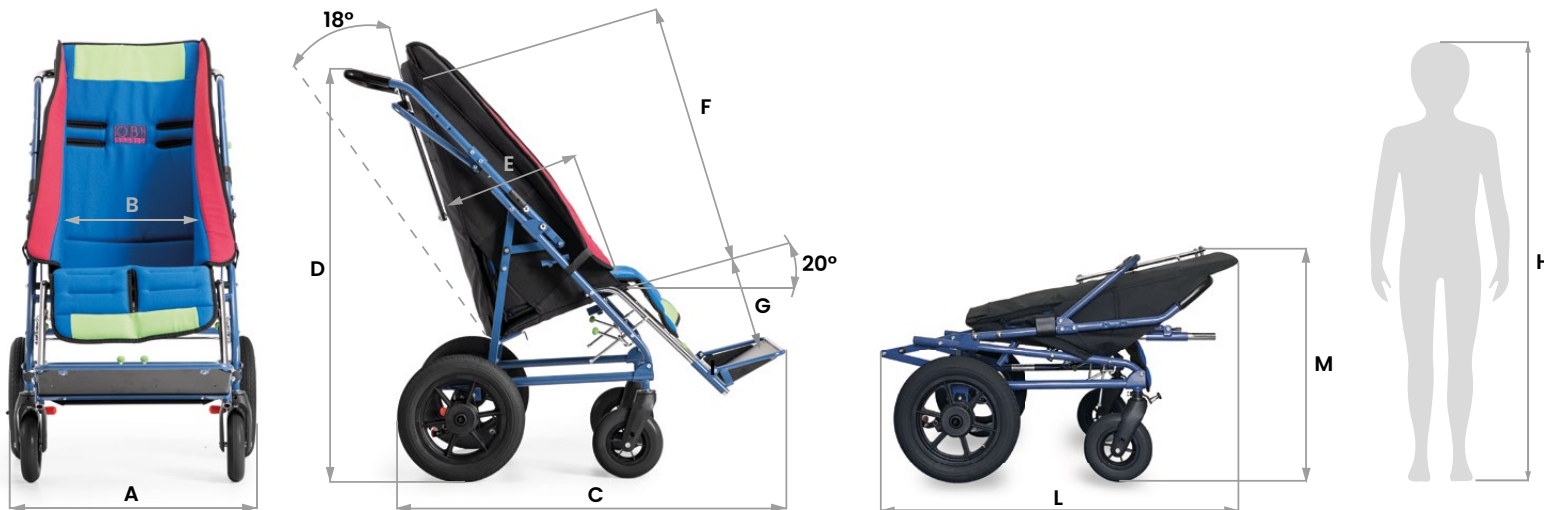





ORMESA[®]
• MADE IN ITALY SINCE 1980 •

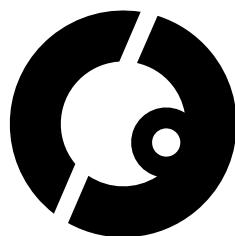

CE complies with MDR (EU) 2017/745

COMPANY WITH
QUALITY SYSTEM
CERTIFIED BY DNV GL
= **ISO 13485** =

Ormesa srl - Via delle Industrie 6/8 | 06034 Foligno (PG) - ITALY
Tel. +39 0742 22927 | FAX +39 0742 22637 | info@ormesa.com

www.ormesa.com