






DISTRIBUTOR/AGENT:	DATE OF ORDER:
	DATE OF DELIVERY:
	NAME OF THE ADDRESSEE/USER:



ACCORDING TO THE ISO  
7176/19 and  
ANSI RESNA WC/19

- size Small      Upholstery colour:  dark grey/green  
 size Medium       dark grey/pink  
 Frame colour: **black**       light grey/green  
     light grey/pink

M  
E  
A  
S  
U  
R  
E  
S

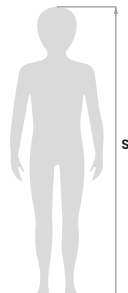	A	B	C	D	E	F	G	J	H	I	L	M	N	O	P	Q	R	S	max. load	frame weight	seat weight	total weight
<b>Small</b>	20-30 cm 7,8-11,8 in	52-62 cm 20,4-24,4 in	37 cm 14,5 in	52 cm 20,4 in	20-30 cm 7,8-11,8 in	14-30 cm 5,5-11,8 in	0°-35°	35°	100-110 cm 39,3-43,3 in	57 cm 22,4 in	67-96 cm 26,3-37,7 in	35 cm 13,7 in	88 cm 34,6 in	33 cm 12,9 in	82 cm 32,2 in	50 cm 19,6 in	81 cm 31,8 in	90-120 cm 35,4-47,2 in	35 kg 77,1 lb	9,5 kg 20,9 lb	8,5 kg 18,7 lb	18 kg 39,6 in
<b>Medium</b>	25-36 cm 9,8-14,1 in	60-75 cm 23,6-29,5 in	46 cm 18,1 in	55 cm 21,6 in	30-40 cm 11,8-15,7 in	20-38 cm 7,8-14,9 in	5°-30°	40°	112-127 cm 44-50 in	62 cm 24,4 in	80-107 cm 31,4-42,1 in	37 cm 14,5 in	112 cm 44 in	39 cm 15,3 in	94 cm 37 in	52 cm 20,4 in	94 cm 37 in	115-140 cm 45,2-55,1 in	45 kg 99,2 lb	11 kg 24,2 lb	12,7 kg 27,9 lb	23,7 kg 52,2 lb

Hyphen-separated numbers indicate the minimum and maximum range of adjustment

### SEAT UNIT, the only one with a tilt-in-space system incorporated in the seat


- HEADREST incorporated in the backrest: your child's head is always supported during the tilt adjustments of the seat and of the backrest, also without accessories.
- Rigid BACKREST, reclinable and adjustable in height. The hinge coincides with the child's pelvis physiological fulcrum.
- Tilt-in-space SEAT, adjustable in width and depth. Reversible (size Small, base 869).
- Fireproof, recyclable FABRIC. Colours durable and highly resistant to sweat, light, dry cleaning and washing. Maximum resistance to rubbing. Oeko-Tex Certificate - for textile products intended for all fabrics that are in direct contact with the skin, a guarantee of quality and safety for health.
- Tip-up FOOTREST, flexion-extension adjustable, with legrest adjustable in height and tilt; the hinge coincides with the child's knee physiological fulcrum.

### 869 OUTDOOR four wheels base - CRASH TESTED ISO 7176/19 and ANSI RESNA WC/19

- Foldable aluminium FRAME. Push-handle adjustable in tilt: it adapts to the mother's height and reduces the encumbrance of the pushchair.
- Puncture proof WHEELS, light and easy to push on every surface: front 17 cm with directional locks; rear 25 cm (size Small) or 30 cm (size Medium), with shock absorbers and independent drum brakes.

### 856 - 851 indoor high-low bases, 850 indoor compact floor base

- Foldable FRAME in steel. Adjustable in height. It is easy and safe to use, thanks to a gas spring mechanism, hand lever (base 851) and foot pedal (base 856)
- 10 cm swivelling WHEELS with brakes (for base 856)
- 7.5 cm swivelling WHEELS with brakes (for base 851)



**Attention!**  
Read carefully the user and  
maintenance manual

ADDITIONAL COMPONENTS

**863**  
Multiadjustable  
Headrest with  
parietal supports



**852**  
Headrest with  
parietal supports



**942**  
Ergonomic headrest  
for Bug Small



**942**  
Ergonomic headrest  
for Bug Medium



**835**  
Shaped headrest



**934**  
Suspension  
headrest  
attachment. Only  
available for 863  
headrest



**838**  
Adjustable lateral  
supports



**868**  
Wrappable and  
flexible trunk  
supports



**834**  
Padded  
abduction block



**834R**  
Adjustable  
padded abduction  
block



**853**  
Vest harness



**903**  
Five point vest  
harness



**906**  
Five point harness



**894**  
45° pelvic belt




**920**  
4 point belt




**827**  
Foot straps



**839**  
Adjustable front  
handle

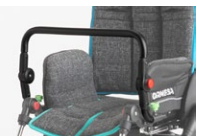

**824**  
Adjustable  
transparent tray

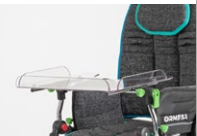

**824B**  
Wide tray only for  
850-851-856 indoor bases.



**905**  
Hand brake  
for 869  
base



**935**  
Oxygen bottle  
holder

Maximum dimensions of the  
oxygen bottle: height 600 mm  
(valve included) – diameter 140  
mm.



**819** Canopy  
 Black/green  Black/pink



**825**  
Rain cover (to fix at the canopy  
819) for 869 base



**818**  
Thermal cover



**858**  
Net basket for 869 base



**911**  
Vent tray



**892**  
Padded covering  
for footrest



**891**  
Set of 4 tie down hooks for the  
forward-facing transport of the  
pushchair in a motor vehicle  
(private/cars, buses...). Only for 869 base



CONFORME ISO 7176/19 e  
ANSI RESNA WC/19



**922**  
Bottle holder for 869 base



**851**  
Indoor high-low  
compact base



**856**  
Indoor high-low  
feeding base

