

The alternative to being bedridden



Marked **CE**

Tie down hooks (4 rings) for the front facing transport of the user in moving vehicles (private cars, conformed buses etc.)  
**AS STANDARD**



DISTRIBUTOR/AGENT:	DATE OF ORDER:
	DATE OF DELIVERY:
	NAME OF THE ADDRESSEE/USER:

**FEATURES**

- **FRAME** made of chromed steel, painted with non-toxic epoxy powder paint. Narrow base to facilitate entering elevators.
- **The seat and backrest are covered with visco-elastic fire-proof material. It is resistant, breathable and anti-slip. It is made of a special and exclusive ORMESA formula:**
  - 1) It adapts to the user's body shapes and it prevents its sliding, with the result of an **EXCELLENT COMFORT and KEEPING POSTURE**, even when changing the backrest recline and tilt (it adapts to the user's body shapes maintaining its support function).
  - 2) Easily washable to keep the **HYGIENE** of the user; it lets **AIR PASS** keeping the **OPTIMAL MICROCLIMATE that HELPS PREVENTING BED SORES**.
  - 3) Durable and fireproof, **GUARANTEE OF DURABILITY and SAFETY**.
  - 4) Ergonomic, breathable, non-toxic, fireproof **UPHOLSTERY**; It does not interfere with the benefits of the visco-elastic cover. The upholstery is fixed with buttons, easy to remove for machine washing.
- Ergonomic **BACKREST**, it holds the shapes of the user. **The special structure of the backrest gives a proper alignment to the trunk and it doesn't generate any pressure on the spine.** Adjustable in tilt by levers with safety pin. The pin allows blocking the levers to maintain the proper adjustment for the user, height adjustable, folding for transportation.
- **Low SEAT** and slightly inclined. Adjustable in depth on four positions.
- **Wide and wrap HEADREST**, padded, multi adjustable (in height, tilt, depth, lateral traslation); thus assuring the **proper alignment between head and trunk**, without interfering with the user's view. Anti-bacterial and Hygienic.
- **ARMRESTS**, padded, easily and stably adjustable in height thanks to the lever.
- **PUSH HANDLE** adjustable in tilt/ height; it helps the care.
- **CONTINUOUS AND WITHOUT EFFORT ADJUSTMENTS**: tilting, backrest recline and legrest tilt are operated by independent gas springs.
- **SIDES** with scissor opening. Thanks to the special joint, the armrests are positioned behind the backrest, eliminating completely all obstacles to the transfer.
- **LEG RESTS**: ergonomic, removable, turnable and adjustable in tilt without effort, activated with gas springs.
- **CALFRETS** adjustable in depth and free tilt to accommodate the movements of the leg.
- **FOOT RESTS**, adjustable in flexion-extension and height.
- **8" front solid swiveling wheels and 12" rear wheels**; drum brakes are operated by a single control pedal.

Please check off the interested products:

\_\_\_\_\_


**COLORS MEASURES**

FRAME COLOUR:  white  black

UPHOLSTERY COLOUR:  dark grey  light grey

SIZES:  36  40  45  50

**Attention!**  
Read carefully the user and maintenance manual

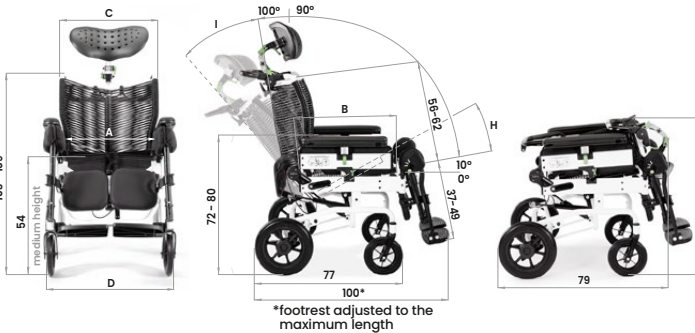

Juditta B30	size 36	size 40	size 45	size 50
<b>A seat bottom width</b>	cm 36	cm 40	cm 45	cm 50
<b>B seat bottom depth</b>	cm 40-43-46-49	cm 40-43-46-49	cm 40-43-46-49	cm 40-43-46-49
<b>C backrest width</b>	cm 47	cm 47	cm 52	cm 57
<b>D overall width</b>	cm 62	cm 62	cm 67	cm 72
<b>Hseat bottom inclination-tilt</b>	from 10° to 35°	from 10° to 35°	from 10° to 35°	from 10° to 35°
<b>I backrest inclination</b>	from 10° to 38°	from 10° to 38°	from 10° to 38°	from 10° to 38°
<b>Total weight</b>	kg 38	kg 38	kg 41	kg 42
<b>Max capacity</b>	kg 120	kg 120	kg 130	kg 140

Hyphen separated numbers indicate a minimum and maximum range of adjustment

**ADDITIONAL COMPONENTS**

<input type="checkbox"/> <b>824-R</b> Reclinable tray, adjustable in tilt, backwards-forwards, with security block  It is made by resistant and hygienic material, design with no edges, it doesn't keep residual food.  	<input type="checkbox"/> <b>894 45°</b> Pelvic belt  	<input type="checkbox"/> <b>834 Padded</b> abduction block  	<input type="checkbox"/> <b>892</b> Footplate covering  	<input type="checkbox"/> <b>905</b> Set of hand brakes with single lever  
<input type="checkbox"/> <b>914</b> Footrest lock  	<input type="checkbox"/> <b>896</b> Tray covering  	<input type="checkbox"/> <b>916</b> Padded leg panel  	<input type="checkbox"/> <b>932</b> Anti-tip device  	<input type="checkbox"/> <b>933</b> IV Pole  
		<input type="checkbox"/> <b>936</b> Headrest cover  	<input type="checkbox"/> <b>938</b> Water-proof cover  	